NAUSICAA

MEDICAL

www.nausicaa-medical.com

WAY UP'GO

FRANCAIS

ENGLISH

FR

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MAXIMUM CAPACITY: 180 kg

CLASS I MEDICAL DEVICE

ALL OUR DEVICES COMPLY WITH THE REQUIREMENTS OF STANDARD NF EN ISO :

10535 : 2021

Contents

Assembly guide / Stand up lift: WAY UP'GO

• Presentation, contraindications, contents	33
• Assembly instructions	34-35
• Using the stand up lift	36-39
• Using the motor	40-42
• Using the strap	43-44
• Storing the lift	45-46
• Re-use of the lift / labelling	47
• Labelling	47-49
• General safety instructions	50
• Cleaning and care	51
• Preventive maintenance and safety checks	52-53
• Troubleshooting guide	54
• Technical specifications and dimensions	55
• Spare parts	56
• Motor technical specifications	57-58
• Warranty	59

Way Up'Go

NAUSICAA
MEDICAL

The WAY UP GO stand up lift is a Class I medical device in accordance with Regulation (EU) 2017/745.

Presentation

The WAY UP GO is a device designed for use in hospitals, clinics, retirement homes, nursing homes, at home... and must be operated by a caregiver. It is designed for use with an appropriate sling. This device has been designed for patients weighing up to 180 kg.

The climatic conditions for using the WAY UP GO are as follows: ambient temperature from 0°C to 40°C, humidity from 20% to 80% and air pressure from 700 to 1060 hPa.

The device can be used in the bathroom (without contact with water) or the toilet. Standing lifts enable the patient to be actively transferred from a bed, a recliner, a toilet or a wheelchair.

Contraindications

Conditions such as osteogenesis imperfecta, osteoporosis, spinal column disease, mental illness or epilepsy may be contraindications.

Contents

The items listed in the following table are included with the lift:

Cardboard	Description	Quantity
WAY UP'GO COMPACT	Stand up lift	1
	Control unit with removable, non-integrated battery	1
	Remote control	1
	User manual	1
	Knee support and knee strap	1
	External battery charger	1

Assembly instructions


THE UNIT IS SUPPLIED FOLDED IN HALF WITH THE KNEE SUPPORT REMOVED.

Step 1 : Loosen the top knobs to release the mast without removing it.

Step 2 : Fully unscrew the lower knobs and remove the mast axle.

Step 3 : Lift the upper part of the unit by twisting it.

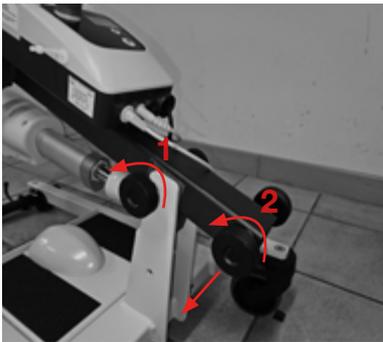
Assembly instructions

Step 4 : Replace the axle, washers and knurls in the lower hole.

Step 5 : Tighten the four knurls to secure the assembly.


Step 6 : Insert the knee support into its housing and tighten its knob to secure.

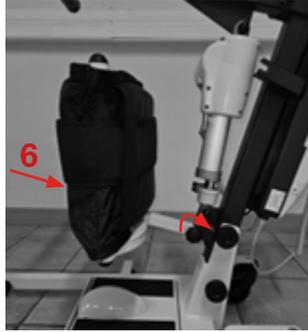Using the device

BEFORE USING THE LIFT, CHECK THAT :

- The wheels turn and roll correctly.
- The brakes work properly.
- There is no wear or deformation on the hooks.
- The various components are correctly connected.
- The battery is sufficiently charged.

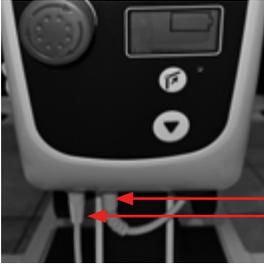Cylinder connector
Remote control connector

Motorised wheel connector

Battery level display

LCD display

Using the device

- Your stand up lift is designed for transferring a patient and must not be used for any other purpose or for moving a patient over long distances.
- Check that the patient's weight does not exceed the maximum weight the lift can support.
- Manoeuvre the stand up lift by pushing the handlebars, never by pushing the patient.
- When transferring a patient, the lift must be handled with care and at a speed appropriate to the situation.
- Use the lift on level, smooth surfaces. It is not recommended to use the lift on a slope of more than 5°. If you have to use a ramp, it is advisable to have a second person to help you. Never use a stand up lift under a shower.

Never recharge the batteries of a stand up lift near a bath or shower.

- Stand up lift are useful and effective when a patient has a certain degree of dependence and needs a caregiver to carry out basic care.
- They help to stimulate the patient and improve mobility; stimulate the cardiac system; combat osteoporosis; and improve the quality of life; Osteoporosis and problems associated with immobility, such as falls during transfers; Stimulate brain activity, helping to maintain urinary and faecal continence; improve intestinal and bowel activity ; Intestinal activity and bladder function.


CAUTION :
DO NOT USE THE DEVICE WITHOUT
THE KNEE SUPPORT.

When using a stand up lift, patients can be divided into two categories:

1. Patients with low muscle tone
2. Patients with high muscle tone

Using the device

- Patients with **low muscle tone** should be transferred to a **semi-seated position**.

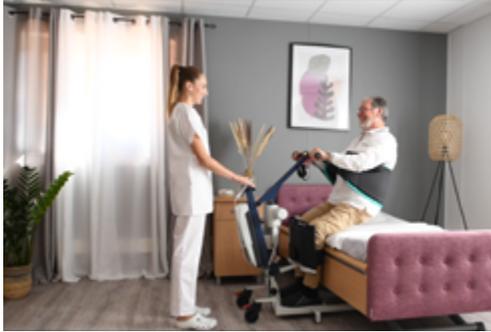- Patients with **medium muscle tone** should be transferred to a **standing position**

- **For safe and effective raising, the patient's hands should be placed on the handlebars on the handlebars provided (see photos above).**

- The patient's shins must be in contact with the knee support at all times.
- The knee support should be positioned by a caregiver.
- Depending on the size of the patient, adjust the position of the knee support strap by sliding it on its support.

Using the device

Braking system :

The brakes are one of the important safety features of the lift. The rear wheels are equipped with brakes that act on the wheels. To lock the wheels, press down on the red part of the wheel. until it stops, using your foot (braked wheel). To unlock the wheels, push the red part upwards

When locking the wheels, always ensure that both wheels are locked.

If only one wheel is locked on a slope, the unlocked wheel will roll and the lift may tip over.

If only one wheel is locked on a slope, the unlocked wheel will roll away and the lift may tip over on its side.

Double wheel with brake on tight.


Double wheel with brake released

Using the motorisation

The WAY UP GO COMPACT has two separate drives:

- An electric actuator to raise the patient.
- A motorised wheel to facilitate patient transfer.

These two functions are controlled by a lithium battery-operated control box:

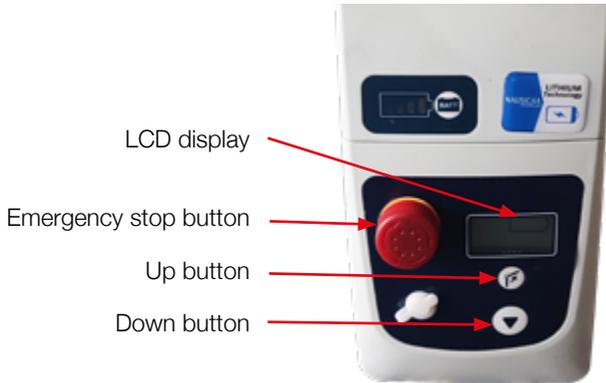Using the motorisation

Switching the control unit on and off :

- To switch the control unit on, turn the emergency stop button a quarter turn. The LCD screen appears, indicating the battery charge level.
- Press the emergency stop button after use or in the event of imminent danger.

Up/down :

- Press and hold the up or down button to activate the arm. The arm stops automatically. The arm stops automatically when the up or down position is reached.

Motorised travel assistance :

- The WAY UP GO is equipped with a motorised wheel to assist the user during patient ride only.
- This function can only be used when a patient is positioned with their feet on the lift. Without this condition, the motorised wheel rotates in a vacuum and cannot be used.
- In addition, the remote control must be positioned on its support in order to ensure that the caregiver can hold the lift throughout the transfer. The forward/backward buttons are ideally placed under the thumb of the carer's right hand to facilitate easier handling. If the remote control is not correctly positioned on its holder, the function will be deactivated and unusable.

Remote control holder

Remote control properly positioned

Using the motorization

Charging the battery :

Battery is loaded directly via the control box. To connect the external charger supplied with your lift to the socket on the front of the control box. When fully charged, close the plug and store the battery charger in the back pocket of the knee support.

Removing/replacing the battery :

Wall chargers are available as optional equipment.

To remove the battery, press the grey button and remove the battery pack by pulling upwards.


To reinstall the battery, press the grey button and position the battery pack so that the bottom of the battery fits into its housing and the top hook holds it in place.

Using the slings

Nausicaa Médical recommends using our PREMIUM and VERTICONFORT slings. For advice on choosing the right product for your needs, please contact our sales department. To select a sling size, the caregiver must match the chest circumference (sub-axillary sling) or waist circumference (back sling) with the dimensions on the chart below.

SIZE	S	M	L	XL
CHEST SIZE	70/90	90/115	115/140	140/165
BETWEEN WAIST AND HIPS	70/90	90/115	115/140	140/165

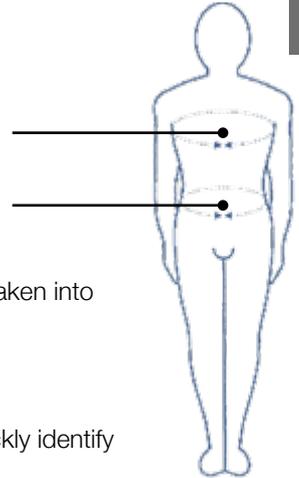ENGLISH

Before using a sit-to-stand lift, the following parameters should be taken into account : the physical handicap, the patient's pathology and general morphology.

Our R&D department has introduced a colour code to help you quickly identify the size of each strap:

- size S is yellow
- size M is red
- size L is green
- size XL is blue
- size TU (one size fits all) is black

TU	A	105
	B	19

IMPORTANT RECOMMENDATIONS

- TO ENSURE MAXIMUM EFFECTIVENESS WHEN USING THE PRODUCTS IN THIS RANGE, IT IS ESSENTIAL TO :
- CHOOSE THE RIGHT SIZE FOR THE PATIENT
- FIT THE PRODUCT AS SNUGLY AS POSSIBLE ON THE PATIENT
- THESE PRODUCTS MUST NOT BE PLACED IN DIRECT CONTACT WITH DAMAGED SKIN.
- READ THE INSTRUCTIONS FOR USE BEFORE USING THE SLINGS.
- BEFORE TRANSFERRING A PATIENT WITH THE STAND UP LIFT, ENSURE THAT THE SLING IS SECURELY HOOKED TO BOTH HOOKS ON THE LIFT.

Using the sling

Position the bottom of the backrest at lumbar level so that the the sling hooks are accessible from either side of the patient.

Pass the backrest straps under the patient's arms.

Attach the two straps to the hooks on the stand up lift.

Adjustments

The multiple-hanger system allows the patient's position to be adjusted.

Vertical positioning of the patient:

The sling is positioned at waist level, and the hooking adjustment is made with the sling under tension and the patient's arms positioned à 120°.

The hands are placed on the handles and the arms are stretched out for more chest tonicity

A quick and full extension of the lift arm produces a standing position.

Partial deployment of the lift arm combined with the under-boot sling gives a seated position.

These slings help to secure the patient while walking.

Reminder of standard ISO 10535 : 2021 on inspection :

B.2.4 Inspection of body support elements (= slings)

A periodic check of the slings should be carried out at the frequency indicated by the manufacturer, i.e. at least twice a year.

More frequent checks may be necessary if a sling is used or cleaned more frequently than the manufacturer's instructions.

Inspections should be carried out by a qualified person who is familiar with the maintenance of the slings.

The inspection should enable to detect signs of deterioration, wear or potential failure.

The inspection report should be kept in a secure place so that it can be consulted in the event of an incident.

This record should include the following information:

- the date of the inspection
- the identification details of the device and its serial number;
- information on the condition of the device;
- the next scheduled inspection date;
- the name, contact details and signature of the inspector.

Storage


To optimise storage, the WAY UP GO COMPACT is designed to be folded for storage. To do this, follow the steps described below :

Step 1: Raise the arm to the top position using **the remote control**.

Step 2: Unscrew the knee support and insert it into the top bracket. Tighten the knob after **inserting the knee support**.

Step 3: Loosen the knobs at the top of the base without removing them.


Unscrew the knobs at the bottom of the base and remove the axle **while holding the upper part of the base in place**.

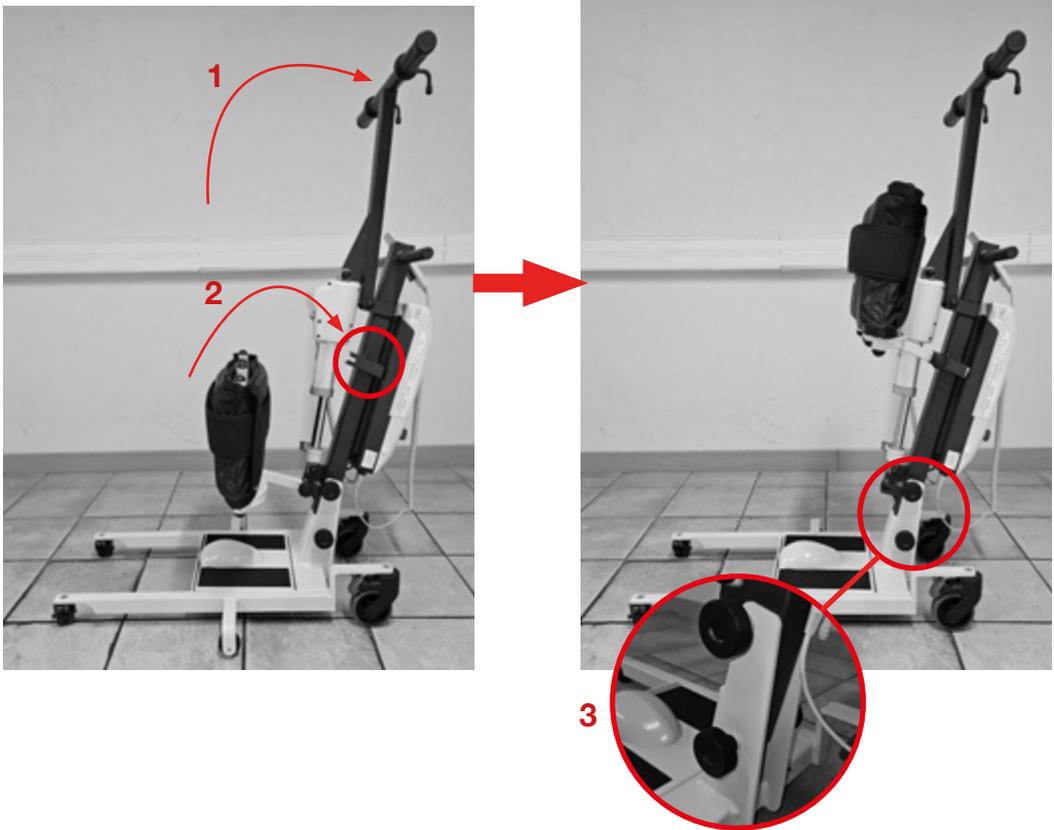Storage

Step 4 : Tilt the top section forward to fold the lift in half, then screw the pin then screw the axle, washers and knurls back onto the base of the mast.

Step 5 : Place the knee support strap under the frame of the device to keep it folded.

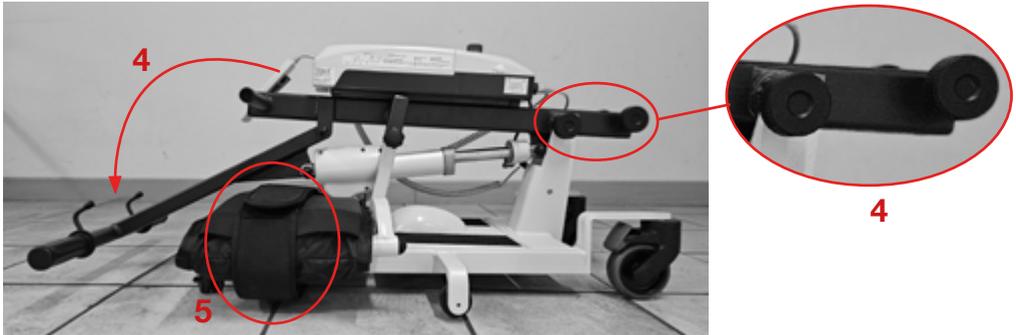- The lift can now be raised using the arm as a handle.
- It can be rolled like a trolley or caddy to take it to its storage location.
- The unit can be stored in an upright position (an optional cover is available).

Dimensions L*W*H = 590*490*1100

Re-use of the stand up lift

The lift can be re-used.


The lift and slings should be cleaned, disinfected and maintained after each use.

Before any re-use, please refer to the chapters :

- Preventive Maintenance & Safety Checks
- Cleaning & Care

Labelling

Labelling on the carton

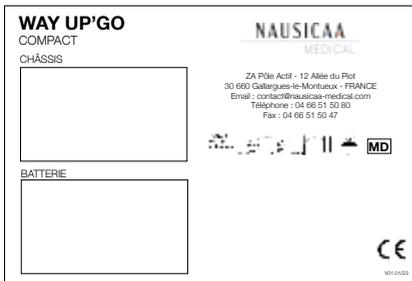Labelling

1

Way Up'Go

NAUSICAA
MEDICAL

2

 SAFETY INSTRUCTIONS - Before each use

NAUSICAA
MEDICAL

- Check the perfect condition of the hoist sling.
- The sling must never be used to move the hoist.
- Check the perfect condition of the traction sling (Easylev range).
- Verify the presence and conditions of safety devices : all axis and their constrained fixations (pins, clips and bolts) as well as the tightness of the screws (see user manuals).

DO NOT USE THE DEVICE IN CASE OF ANY LACK OR BAD OVERALL STATE

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 BATTERY MAINTENANCE

NAUSICAA
MEDICAL

A POORLY MAINTAINED BATTERY DETERIORATES

- Regularly charge the battery at the maximum.
- Batteries must be used and stored at temperatures between 5°C and 40°C.
- The charger should be checked if the battery stops working.
- Making small incomplete charges shortens the battery life.

LEAD BATTERY

- Full charge twice a week when used at home.
- Battery lifetime depends on its correct loading.

LITHIUM BATTERY

- Can be recharged as often as desired.
- If the charge level indicator remains off, recharge the battery.

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Labelling

General safety instructions

1. Use the lift only for its intended purpose, in compliance with current legislation on medical devices, the relevant legal regulations, the regulations for the protection of workers and accident prevention regulations.
2. Note that the lift is a medical device and that the user must therefore comply with the regulations governing the use of medical devices.
3. The electrical installation requirements for the room or area where the stand up lift is used must be controlled.
- 4) The lift should only be used after the patient has been trained in its use and with full knowledge of the device..
5. Before using the lift, read all instructions in order to avoid damage caused by improper handling or exposure to hazards. The instruction guide contains important information and notes that are necessary for the use of the lift and can be accessed via the QR code on the lift.
6. Before use, check that the lift and its accessories are in good conditions.
7. Before using the stander with other medical or non-medical devices, check that the combination of these products is authorized and that they can be used together safely.
8. The assembly, commissioning, maintenance and repair of the stander shall be entrusted only to the
to competent specialists.
9. It is the responsibility of the user or operator to ensure (by means of appropriate measures and instructions) that a Mechanical stress of the of the load cord (by bending, pulling, shearing, crushing) is excluded
while loading or cleaning the part. This also applies to the electrical cables of others
devices used with the stander.
10. Observe the start-up time and the maximum permissible load. These values should not be outdated, otherwise safe operation will no longer be ensured.
11. Do not expose the Stander to direct sunlight, neither to heat nor to humidity.
12. Make sure that no moisture enters the electrical system.
13. Avoid mechanical stress on the electrical cords used. Pulling, bending, or crushing cords may damage them.
14. Recharge batteries in a well-ventilated area
15. Electromagnetic or other disturbances cannot be excluded between the stander and other Devices. If there is a risk of such disturbances, the source of the disturbance must be far away or the sitter should not be used.
16. Disruptions caused by the use of mobile communication devices cannot be completely excluded. For this reason, a safety distance of at least 3.3 m must be maintained in order to ensure the safety of the works of at least 3.3 m in order to ensure the
Safe operation of the stander.
17. Do not allow children to stand unsupervised near the stander.
18. The stander should no longer be used as soon as abnormal noises, damage or any other another malfunction. In this case, do not connect the stander to the charger, but inform NAUSICAA MEDICAL
19. If damaged or defective, the stander should no longer be used and should no longer be connected to the
sector. Inform the dealer to instruct him to remedy the defect or failure.
20. Any serious incident that occurs in connection with the device should be reported to NAUSICAA
MEDICAL and the competent authority of the Member State in which the user and/or patient is established

Cleaning and maintenance

BEFORE ANY MANIPULATION:

- UNPLUG THE POWER CHARGER FROM THE MAINS.
- CHECK THAT ALL ELECTRICAL ELEMENTS ARE CONNECTED TO EACH OTHER.

DAILY MAINTENANCE / BETWEEN EACH PATIENT:

- Clean with a surface product.
- Clean the cylinder casing and cables, if they have been splashed with fluids bodily infections, including urine.
- Clean and disinfect all parts in direct contact with the patient (handles, foam)

COMPLETE CLEANING AND DISINFECTION:

- Clean the entire chassis of the appliance with a pre-disinfectant detergent.
- Clean the engine with a cloth slightly dampened with a pre-disinfectant detergent solution
- Never clean the appliance, especially the electrical system, with a high-pressure cleaner, a water hose or similar.
- Disinfect the appliance (by spraying with surface disinfectant or by misting of disinfectant solution for example), please respect the contact times of the product used.

MAINTENANCE AFTER CLEANING AND DISINFECTION :

- Check the appliance for visible damage.
- Check that no parts are missing.
- Check that the wheels are working properly and that there is no material staining their bearing (hair, scraps of carpet, etc.).
- Check that the controls are working properly as well as that the remote control is connected, and of the cylinder on the battery pack.
- Check the integrity of the electrical cables (cylinder and remote control)

CAREFUL:

- **THE DETERGENTS USED MUST HAVE A NEUTRAL PH.**
- **AVOID ABRASIVE PRODUCTS AND SOLVENT, THEY COULD DAMAGE THE SURFACE FINISH OF THE DEVICE.**

If you notice a fault, contact the after sales service:

Stand Up Lifts & Hoists
Phone : 04 66 71 71 80
Fax: 04 66 71 71 81
Mail: sav@nausicaa-medical.com

Preventive maintenance and safety checks

Stand Up type: WAY UP'GO COMPACT

Chassis serial number :

Maintenance must be carried out by trained personnel. Checking safety points
(frequency: depending on use, at least once a year) :

This check is carried out visually; any sign of damage must be replaced.

ENGLISH

		COMPLIANT	NOT COMPLIANT	GREASING	DATE OF CHANGE
1	Cylinder upper clevis mounting (bolted)				
2	Lower cylinder clevis mounting (bolted)				
3	Mast / lifting arm assembly fastening (bolted)				
4	Wheel pivot attachment (bolted)				
5	Fixing the base plate to the mast (knurled)			X	
6	Knee rest attachment (knurled)			X	
7	Structural condition (welds, joints, oxidation)			X	
8	Fixing the motorised wheel			X	
9	Condition of cables and connectors			X	
10	Check that labels are present and legible			X	
11	Check the platform clearance and its knocks			X	

*use a 3-in-1 product or similar



Electrical functions check (frequency: according to use, at least once a year):

This check is carried out with the stand up lift under load.

		COMPLIANT	NON COMPLIANT	DATE OF CHANGE
A	Lifting cylinder operation (up/down)			
B	Emergency stop operation			
C	Battery status (charge conservation)			
D	Remote control status			
E	Motorised wheel operation (forward/reverse)			
F	Brake operation			
G	Absence of suspicious noise			

Checks carried out on :	Checks carried out by :	Date of next inspection

Preventive maintenance and safety checks

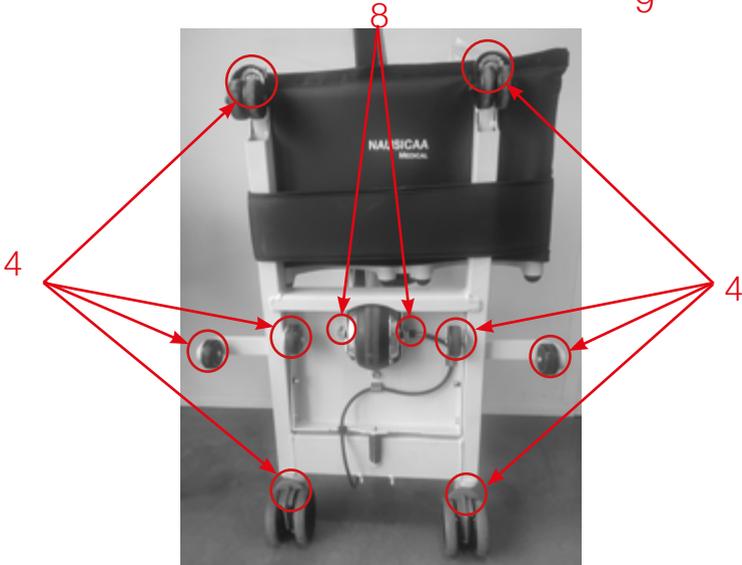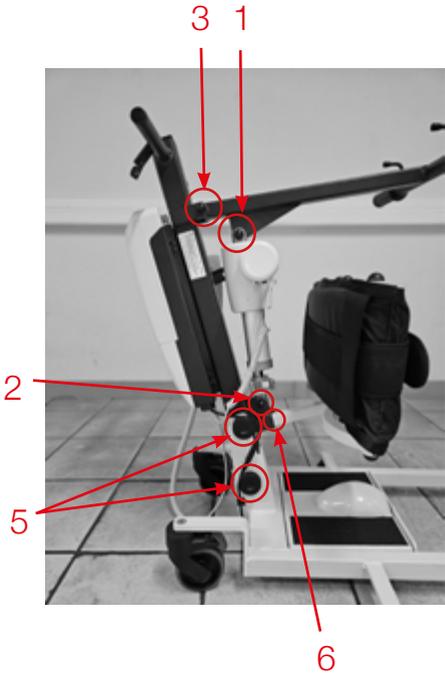Troubleshooting guide

Common faults :

Symptoms	Causes	Solutions
The appliance squeaks and moving parts	Lack of lubrication	Lubricate moving parts
The jack does not work but you can hear a click in the control box when the remote control is operated	<ol style="list-style-type: none"> 1- The battery is discharged 2- The jack cable is disconnected or incorrectly connected 3- The jack or the control box are faulty 	<ol style="list-style-type: none"> 1- Recharge the battery 2- Reconnect the jack cable 3- Have the electrical assembly serviced
The jack does not work and there is no click in the control box	<ol style="list-style-type: none"> 1- The remote control is not on its bracket 2- The magnet attached to the remote control bracket is missing 3- The motorised wheel cable is disconnected or damaged 4- The control box, remote control or wheel are faulty 5- The emergency stop button is engaged 	<ol style="list-style-type: none"> 1- Position the remote control correctly on its holder before using transfer assistance 2- Have the magnet replaced 3- Check the connection and cable of the motorised wheel 4- Have the electrical system serviced 5- Unlock the emergency stop button
The wheel does not work	<ol style="list-style-type: none"> 1- The emergency stop button is pressed 2- The battery is discharged 3- The battery is not correctly installed 4- The remote control is disconnected 5- Control box or remote control faulty 	<ol style="list-style-type: none"> 1- Unlock the emergency stop button 2- Charge the device 3- Insert the battery correctly 4- Connect the remote control 5- Have the electrical system serviced
The motorised wheel turns but no longer the lift	The mobile platform does not tilt correctly and is blocked	Check that nothing is preventing the platform from tipping over and that the patient's feet are correctly positioned on the platform.

Technical specifications and dimensions


- Steel frames
 - Laser-cut and bent sheet metal, tubes and sections
 - Oven-cured epoxy paint
 - Ball-bearing wheels
 - Total weight of compact version: 32.5 kg
 - Weight of top section: 17.40 kg (with one axle and two wheels)
 - Weight of compact version base: 15.1 kg (with other axle and two wheels)
 - Maximum load: 180 kg
 - Product life: 8 years* (excluding electrical part)
- * Subject to regular maintenance

Dimensional specifications

Dimensions

- A.** Total length: 915 mm
- B.** Minimum useful height: 770 mm
- C.** Maximum useful height: 1160 mm
- D.** Total height: 910 mm
- E.** Frame height: 80 mm
- F.** Minimum base width: 390 mm
- G.** Maximum base width: 590 mm
- H.** Turning diameter: 890 mm
- I.** Ground clearance: 43 mm

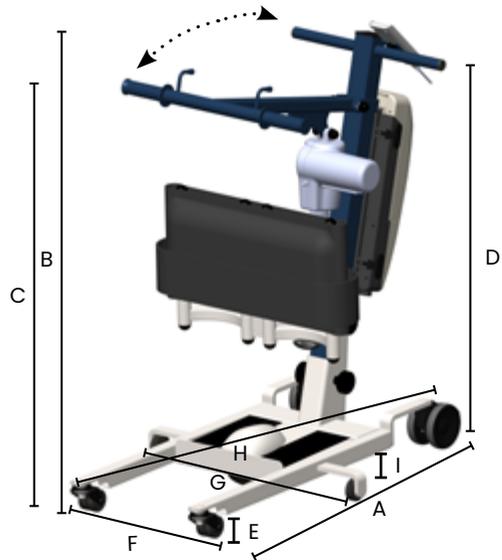Spare Parts

**For all requests for spare parts
please contact our after-sales service:**

Verticals & patient lifts

Telephone: 04 66 71 71 80

Fax: 04 66 71 71 81

Mail: sav@nausicaa-medical.com

Motor specifications

Power supply :

- Removable 24V battery pack (Lithium)
- Digital range display
- Independent wall charger
- IP66 PROTECTION

Class II
equipment

ENGLISH

Electric actuator :

- Electric safety lowering
- Low-voltage DC motor 24V/180VA
- Maximum thrust force 8000N
- 11cm stroke
- Limit stop by contactor
- Noise level: less than 55dB at 1m
- IP55 PROTECTION

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- ABS case
- Electronic overload protection
- Emergency stop button
- Electromagnetic compatibility: complies with standard EN60601-1-2
- Complies with EN60601-1
- Protection class II
- IP54

Remote control :

- Button pressure: less than 1N
- 4 functions
- Spiral cable with low control voltage (5Vdc)
- IP56

Motorised wheel :

- 24V brushless motor
- Maximum torque 5.5N.m
- IP54


Connection scheme

Motorisation technical specifications

Control box display : :

The control box has a display to provide the following information to the user :

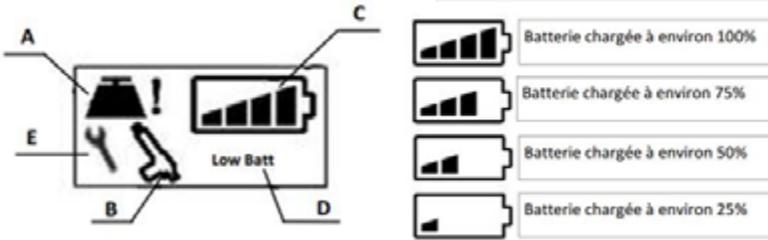A-Overload indicator: The load to be lifted is too high and could damage the device and endanger users. The lift must be relieved.

B-Unnormal operating time indicator: The lift has been used for too long compared to its idle time. This is a safety message to protect the integrity of the lift. The expected operating ratio is 10% use time to 90% rest time. Allow the lift to rest before attempting to use it again.

C-Battery level indicator.

D-Low battery indicator. Risk of damage to the battery, which must be recharged before use. An audible signal warns the user as soon as the battery reaches 25% of its charge.

E-Cylinder maintenance or replacement indicator. The device has been used for a total 600,000 seconds, i.e. approximately 10,000 cycles. A general overhaul of the whole device is necessary.

Guaranty

- Article 1: NAUSICAA Médical S.A.S. guarantees this device against all manufacturing and assembly defects of its mechanical and electrical components, and only for equipment used under the conditions stipulated by NAUSICAA Médical S.A.S.

The warranty covers both mechanical and electrical parts.

This guarantee is valid for 5 years except for the lithium battery which is guaranteed for 2 years following the requirements below.

- Article 2 : Article 2: The warranty entitles the replacement of the defective parts or a replacement by an equivalent device. We do not bear the cost of labour.

- Article 3: The reseller is responsible for the shipping costs of the device and all related costs. The goods always travel at the reseller's risk and responsibility.

Returning costs of the goods to the after sales service will be bear by the distributor or reseller.

- Article 4 : The guarantee does not apply if the complaints are the result of :

- an accident, misuse of the appliance or negligence on the part of the purchaser.
- transport of the appliance without adequate protection.
- modifications or alterations not approved by NAUSICAA Médical S.A.S.
- the impact of external factors (natural disasters, fire, impact, humidity, flooding, lightning, etc.)
- installation and/or use in a manner that does not comply with technical and safety standards, if the device is operated in a country other than the country of purchase; and/or if the power supply is not suitable for the voltage at which the appliance is used.
- failure to carry out routine maintenance.

- Article 5 : The retailer may not invoke the benefit of the guarantee :

- if the serial number of the appliance has been removed, modified or rendered illegible.
- if the equipment under guarantee has been modified without the approval of NAUSICAA Médical S.A.S.

- Article 6 : Article 6: During the repair of defective equipment, no equipment will be loaned.

- Article 7 : All warranty claims must be made through the retailer. If this is not possible, the purchaser may send the equipment directly to NAUSICAA Médical S.A.S.

In this case, the purchaser must indicate, in a letter enclosed with the equipment, the contact details of the retailer and a copy of the certificate of warranty and a copy of the purchase invoice.

- Article 8 : Spare parts under guarantee will only be sent after consultation with the NAUSICAA Médical S.A.S. After Sales Service. It should be noted that defective spare parts must be returned to NAUSICAA Médical S.A.S.

- Article 9 : Defective parts changed under or outside the guarantee will be guaranteed for 6 months from the date of repair or dispatch of spare parts.

- Article 10 : No retailer may unilaterally modify the terms of this guarantee.

NAUSICAA

MEDICAL

EXPORT CONTACT HEAD OFFICE

Telephone : 04 66 51 50 80

Fax : 04 66 51 50 47

Mail : export@nausicaa-medical.com

www.nausicaa-medical.com

Your contact :

manufactured by NAUSICAA Médical S.A.S.